

# **Easy Fit Portable Gantry Data Sheet**

### **Clinical Intention**

The Easy Fit Portable gantry system is a great solution for situations where a permanent ceiling track cannot be installed. It is easy to set up and does not require any mounting on walls, ceilings, or the floor. The gantry is designed to work with a portable lift and sling, and together, they form the Easy Fit portable gantry system. This system makes it possible to move individuals with mobility impairments with minimal strain or risk to the caregiver while ensuring complete safety, dignity and comfort for the person being moved. It is important to note that the gantry is designed for internal use only.

#### **Applied Parts**

The Easy Fit Portable Gantry requires the following components to complete the system.

- **Portable Lift:** The gantry requires a CPP Lift to raise and lower the patient within the gantry track system. For more details, please visit the website. While other portable lifts may be used, it's important to contact your local authorized dealers to obtain approval.
- **Sling:** To view the complete list of applicable slings, please visit the website.



#### **Features & Benefits**

**Easy Assembly:** The Easy Fit portable gantry does not require any tools for installation. It can be assembled and disassembled with ease and is lightweight making it very portable, this allows for simple relocation when required. The gantry can be assembled by a single person.

**Structure Alternation Free:** The gantry is designed to be free-standing; this means that there are no structural alternations required for installation. Huge benefit over a ceiling track system.


**Stay at home:** The gantry can extend the duration individuals can remain in the comfort of their own homes by several years. When combined with a portable lift, it will meet the user's lifting and handling needs with dignity and comfort.

**Fully Adjustable:** The gantry allows for height and width adjustment, allowing it to be set up for your personal requirements. The gantry allows for a portable lift to traverse along the length of its track.

**Mackworth CPP:** The gantry is compatible with our portable Mackworth CPP lift. Together with a sling, they will meet your moving hand handling



## **Technical Specification**

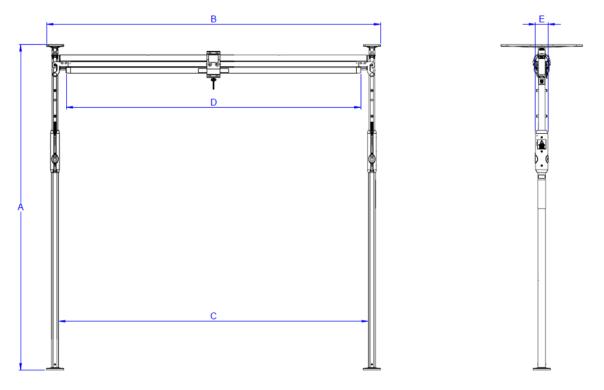

Gantry Dimensions	
A – Minimum - Maximum Height	80.55" – 108.27"
B – Minimum – Maximum Width	82.80" - 117.01"
C – Minimum - Maximum Internal Clearance	76.77" – 111.02"
D – Minimum – Maximum Horizontal Travel	73.07" – 107.32
E – Leg Width	3.15"
General	
User Maximum Weight	440lb
Product Weight	< 50lbs
Expected Lifetime	10 yrs
Construction Type	Aluminium
Packaging Specification	
Total Packaging Weight	< 60lbs
Packaging Box Dimensions (L x W x H)	83.7" x 11.6" x 5.51"
Applied Standards	
IEC 62366-1:2015:	BS EN ISO 14971:2019
BS EN ISO 10535:2006	
Certifications & Accreditations	
	CE



Manufactured by: Prism Healthcare (UK)

